

Dry Bath Incubator

GENAXY MK-10/MK-20

The Dry Bath Incubator is a microprocessor-controlled product using advanced thermoelectric technique. With this technique, more stable and accurate temperature control is achieved. Replacement is convenient with optional metal modules of various specifications configured for different types of test tubes.

Features

- ❖ LCD display. It is easy to setup and use
- ❖ Simultaneous display of set and actual time, temperature
- ❖ Temperature deflection adjusting meet with different demands
- ❖ As it is equipped with various optional blocks, the instrument can adapt to different tubes & wells to cope with experimental needs. It is easy to replace the metal blocks and is very simple to clean and sterilize
- ❖ Customized blocks are available to suit your specific demands
- ❖ Microprocessor controlled incubation time and temperature
- ❖ Beep-signal/Stop after program completion
- ❖ Peltier design of MK-20 provides efficient thermal control



GEN-MK-10





GEN-MK-20

Technical Specification

Catalog No	GEN-MK-10	GEN-MK-20
Temperature Control Range	RT.+5°C~100°C	-10°C~100°C
Timing Range	1min~99h59min	1min~99h59min
Temp. Control Accuracy	±0.5°C	±0.5°C
Display Accuracy	±0.1°C	±0.1°C
Heating Time(25 to 100°C)	≤12min	≤15min
Cooling Time(20 to -10°C)	----	≤15min (*)
Heating Part	Heater	Peltier
Cooling Part	----	Peltier
Certificate	CE	CE
Dimension(mm)	300 x 220x 180	300 x 220x 180
Net Weight	2.5kg	5kg

Accessories

Catalog No	Type	Description
AS-01051-01	A	Block, 96 x 0.2ml Tubes
AS-01051-02	B	Block, 54 x 0.5ml Tubes
AS-01051-03	C	Block, 35 x 1.5ml Tubes
AS-01051-04	D	Block, 35 x 2.0ml Tubes
AS-01051-05	E	Block, 15 x 0.5ml + 20 x 1.5ml Tubes
AS-01051-06	F	Block, 24 x Diameter ≤Φ 12mm Tubes
AS-01051-07	G	Block, 12 x 15 ml Falcon Tubes
AS-01051-08	H	Block, 6x50ml Falcon Tubes
AS-01051-09	I	103 x 67x 30 mm Bath Block
AS-01051-10	J	96 x 0.2ml Microplate
AS-01051-11	K	24 x 5ml Tubes
AS-01051-12	L	96-Deep well plate

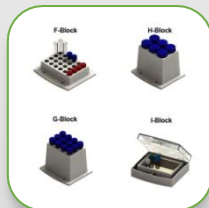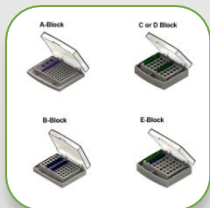

Genaxy Scientific Pvt. Ltd.

155-156, C-Block, Community Centre, Plot No. 32 to 57, Behind Janak Cinema, Janak Puri, New Delhi-110 058

Phones : +91-11-4761 9999 (100 Lines) Fax : +91-11-2553 4163

E-mail : info@genaxy.com