

Gel Casters for SE600 and SE400 Series Vertical Electrophoresis Units

Enables several gradient or homogenous gels to be cast simultaneously

Multiple Gel Casters

- SE675 casts up to four 16 cm long gels
- SE615 casts up to ten 16 cm long gels
- Easy to use and leak free
- Cast homogenous or gradient gels—Pour homogenous gels directly into the top of the casting chamber; pump gradient gels through an inlet port at the bottom

Ordering Information

Cat. #	Description
SE675	Four Gel Caster Includes Glass Plates, 18 x 16 cm—8 pcs
SE615	Ten Gel Caster Includes Glass Plates, 18 x 16 cm—20 pcs

Includes:

- Caster Body with Gasket
- Face Plate
- Clamps—4 pcs
- Glass Plates
- Filter Sheets—5 pcs
- Wax Paper—100 Sheets, Precut
- Spacer-Mate—Alignment Template
- Thumbscrews—2 pcs
- Filler Plug Set
- Black Vinyl Cap
- Gel Seal
- Bubble Level

Order SE511 combs and sets of spacers separately, see page 30.

Dual Gel Caster

SE6015 Dual Gel Caster is included when you order an SE600 Series unit. It can also be ordered separately.

- Casts two 8, 16, or 24 cm long gels at a time
- Double the capacity with optional divider plates (see page 22)
- Easy to use and leak free
- Glass sandwiches seal into the casting cradle with a simple twist of the cams

Ordering Information

Cat. #	Description
SE6015	Dual Gel Caster


Includes:

- Caster Stand with 2 Cradles and Leveling Feet
- Laminated Gaskets—2 pcs

Order glass plates, SE511 combs and sets of spacers separately, see pages 22 and 30.



SE675

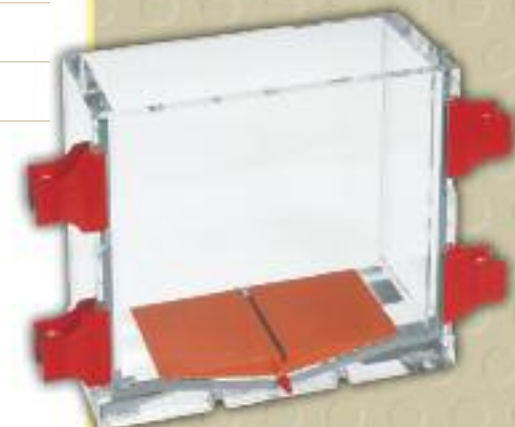

SE615



SE6015