

# Gel Casters for SE600 and SE400 Series Vertical Electrophoresis Units

Enables several gradient or homogenous gels to be cast simultaneously

## Multiple Gel Casters

- SE675 casts up to four 16 cm long gels
- SE615 casts up to ten 16 cm long gels
- Easy to use and leak free
- Cast homogenous or gradient gels—Pour homogenous gels directly into the top of the casting chamber; pump gradient gels through an inlet port at the bottom

## Ordering Information

Cat. #	Description
SE675	Four Gel Caster Includes Glass Plates, 18 x 16 cm—8 pcs
SE615	Ten Gel Caster Includes Glass Plates, 18 x 16 cm—20 pcs

### Includes:

- Caster Body with Gasket
- Face Plate
- Clamps—4 pcs
- Glass Plates
- Filter Sheets—5 pcs
- Wax Paper—100 Sheets, Precut
- Spacer-Mate—Alignment Template
- Thumbscrews—2 pcs
- Filler Plug Set
- Black Vinyl Cap
- Gel Seal
- Bubble Level

Order SE511 combs and sets of spacers separately, see page 30.

## Dual Gel Caster

SE6015 Dual Gel Caster is included when you order an SE600 Series unit. It can also be ordered separately.

- Casts two 8, 16, or 24 cm long gels at a time
- Double the capacity with optional divider plates (see page 22)
- Easy to use and leak free
- Glass sandwiches seal into the casting cradle with a simple twist of the cams

## Ordering Information

Cat. #	Description
SE6015	Dual Gel Caster


### Includes:

- Caster Stand with 2 Cradles and Leveling Feet
- Laminated Gaskets—2 pcs

Order glass plates, SE511 combs and sets of spacers separately, see pages 22 and 30.



SE675

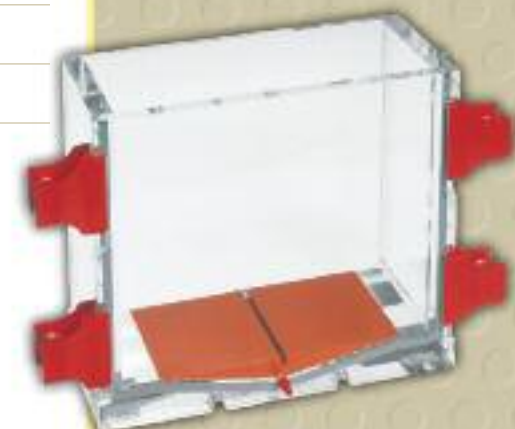

SE615



SE6015